



State of New Hampshire

Department of Safety

Division of Fire Standards and Training & Emergency Medical Services
98 Smokey Bear Boulevard, Concord, New Hampshire
Mailing Address: 33 Hazen Drive, Concord, New Hampshire 03305-0002

Business: (603) 223-4200 Fax: (603) 271-1091



Robert L. Quinn
Commissioner

Justin A. Cutting
Director

CLINICAL BULLETIN

Bulletin #	Title			Date Issued
81	Counterfeited Combat Application Tourniquets			March 1, 2023
Superseded	Released By	Approved By	Source	For More Information
	FSTEMS	Chief Romanello	National Intellectual Property Rights Coordination Center	FST/EMS (603) 223-4200

COUNTERFEIT TOURNIQUET MEDICAL DEVICES

SAFETY HAZARD TO LAW ENFORCEMENT AND FIRST RESPONDERS

Law enforcement has identified a trend regarding the sale of counterfeited Combat Application Tourniquets medical devices. A counterfeit tourniquet is one that bears the trademark of CAT Resources but is produced by another party and is usually not made to the specifications of the original equipment manufacturer. Several public safety agencies have reported failure of the counterfeit tourniquet devices in emergency situations.

ISSUE AND COUNTERFEIT PARTS:

The Combat Application Tourniquet is constructed with a strap that is placed above the wound on the injured appendage, affixed with a Velcro strap, it is tightened into place with a small tension rod (windlass). The tightening of the tension rod applies pressure to restrict blood flow to the affected area. The counterfeit models of the tourniquets are compromised by a poorly constructed tension rod made of plastic or a substandard composite that causes the rod to bend or break, rendering the tourniquet ineffective and potentially deadly as it will not stop the loss of blood.

HOW TO PROTECT YOURSELF:

Only conduct business with the original manufacturer or authorized dealer and use caution when purchasing these medical devices through e-commerce sources. Beware of “too good to be true” prices and use caution when purchasing Combat Application Tourniquets that are shipped from other countries. **Authentic Combat Application Tourniquets are only produced in the United States.**

CHECKING YOUR GEAR:

On pages two and three of this bulletin are photos of genuine and counterfeit tourniquets showing some indicators to look for when checking your gear or purchasing new North American Rescue Combat Application Tourniquets. Check your gear to ensure safety in an emergency medical situation.

WHAT TO DO IF YOU SUSPECT YOU HAVE PURCHASED COUNTERFEITS:

Consult an original equipment manufacturer representative for advice on identifying and replacing the counterfeit medical device. If you believe you have been the victim of a crime, contact your local Homeland Security Investigations (HSI) field office to report the crime or submit an online complaint to the National Intellectual Property Rights Coordination Center at www.iprcenter.gov.



State of New Hampshire Department of Safety

Division of Fire Standards and Training & Emergency Medical Services
98 Smokey Bear Boulevard, Concord, New Hampshire
Mailing Address: 33 Hazen Drive, Concord, New Hampshire 03305-0002


Business: (603) 223-4200 Fax: (603) 271-1091



Robert L. Quinn
Commissioner

Justin A. Cutting
Director

CLINICAL BULLETIN

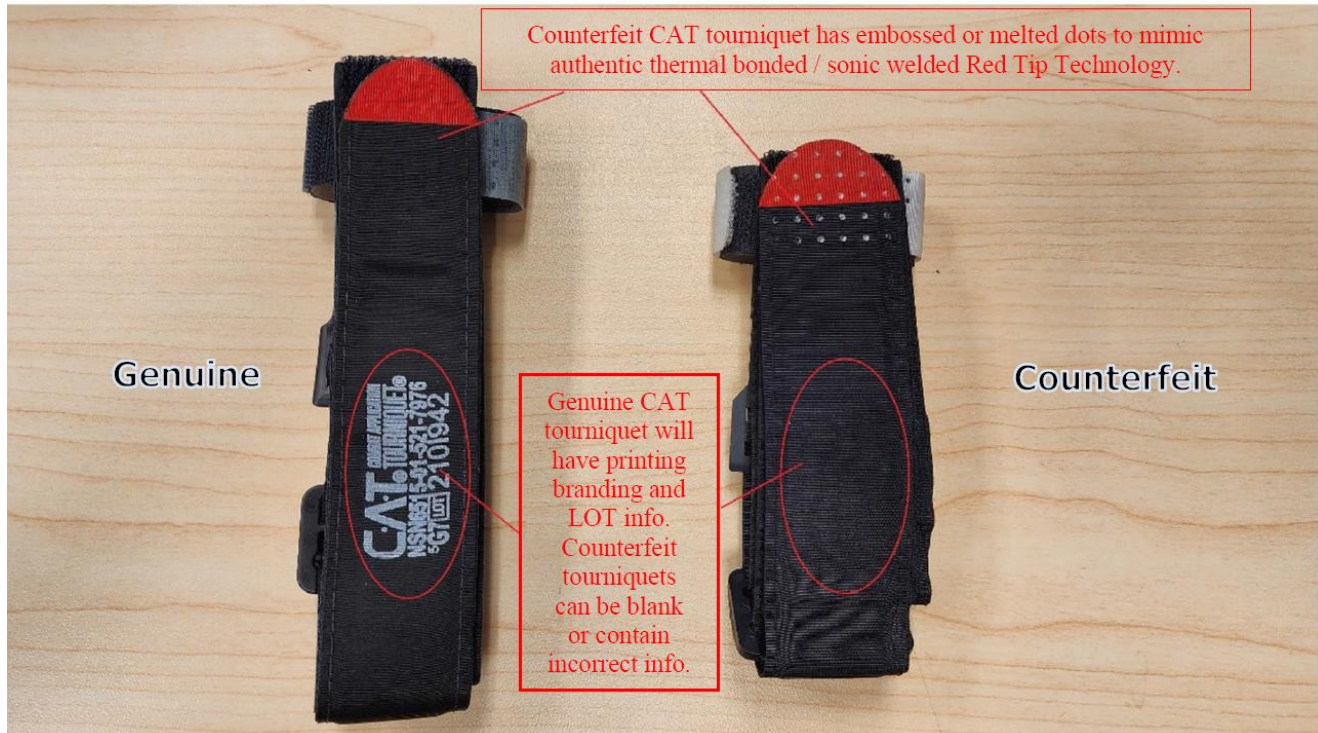



State of New Hampshire

Department of Safety

Division of Fire Standards and Training & Emergency Medical Services
 98 Smokey Bear Boulevard, Concord, New Hampshire
 Mailing Address: 33 Hazen Drive, Concord, New Hampshire 03305-0002

Business: (603) 223-4200 Fax: (603) 271-1091



Robert L. Quinn
 Commissioner

Justin A. Cutting
 Director

CLINICAL BULLETIN

